

Liofilchem®

Certificate of Analysis

Page 1 of 1

Product	Batch	Production date	Expiration date
Ampicillin AMP 0.016-256 µg/mL	020918016	09.02.2018	2021.02.08

Ref.

92003 – 920030 – 920031

Antimicrobial Susceptibility Testing

Tested according to current CLSI methodology

Control strains	Medium	Inoculum	Incubation	Expected Results MIC range (mg/L)	Results MIC(mg/L)
<i>Staphylococcus aureus</i> ATCC® 29213	Mueller Hinton II Agar	0.5 McFarland in saline	35 ± 2°C, ambient 16-20 h	0.5–2	0.5
<i>Enterococcus faecalis</i> ATCC® 29212	Mueller Hinton II Agar	0.5 McFarland in saline	35 ± 2°C, ambient 16-20 h	0.5–2	0.75
<i>Escherichia coli</i> ATCC® 25922	Mueller Hinton II Agar	0.5 McFarland in saline	35 ± 2°C, ambient 16-20 h	2–8	4
<i>Escherichia coli</i> ATCC® 35218	Mueller Hinton II Agar	0.5 McFarland in saline	35 ± 2°C, ambient 16-20 h	>32	>32
<i>Klebsiella pneumoniae</i> ATCC® 700603	Mueller Hinton II Agar	0.5 McFarland in saline	35 ± 2°C, ambient 16-20 h	>128	>256
<i>Haemophilus influenzae</i> ATCC® 49247	Haemophilus Test Agar	0.5 McFarland in broth	35 ± 2°C, 5% CO ₂ 20-24 h	2–8	4
<i>Streptococcus pneumoniae</i> ATCC® 49619	Mueller Hinton II Agar (Sheep blood 5%)	0.5 McFarland in broth	35 ± 2°C, 5% CO ₂ 20-24 h	0.06–0.25	0.064
<i>Haemophilus influenzae</i> ATCC® 49766	Mueller Hinton Fastidious Agar	0.5 McFarland in saline	35 ± 1°C, 5% CO ₂ 16-20 h	0.06–0.25*	0.25

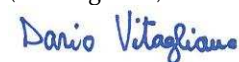
*Established and validated by EUCAST

Batch Release

Approved

Date

26.03.2018

SignatureQuality Control
(D. Vitagliano)

The results reported were obtained at the time of release.